


THE Static Eliminator BY PLAS ■ LABS, INC.

*For Use With
Any Glove Box
Or Isolator!*

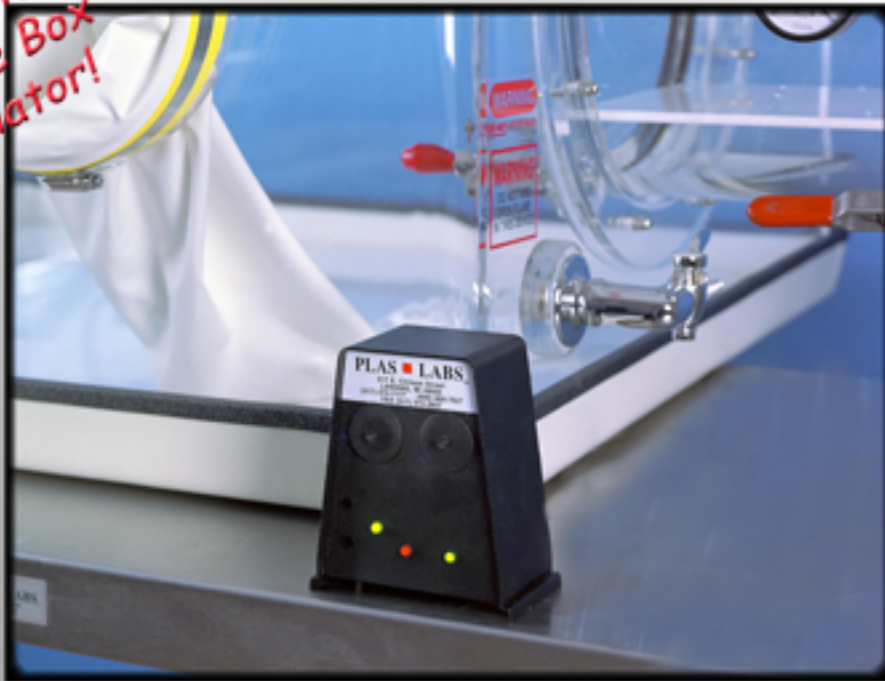

800-AS/SPI

Model #800-AS/SPI, Work Station Ionizer effectively eliminates all static charges within 36" of unit. The operation of the Work Station Ionizer is simple; place the unit on a flat surface and turn the power on. Unit features include: ion pulse rate adjustment, ion balance adjustment and a compact footprint. The Work Station Ionizer is virtually maintenance free; the emitter points should be cleaned when visual accumulation of dust or powders occurs.

Overall Dimensions:	4.25" wide x 2.5" deep x 5.0" high
Input voltage:	110/220 VAC; 50-60 Hz.
Output voltage:	15 kVDC
Electrode (emitter) material:	Thoriated Tungsten (replaceable)
Ion coverage:	up to 36"
Pulse frequency:	Adjustable from 1 pulse/sec. to 5 pulses/sec.
Ozone:	less than 0.1 ppm
Housing:	Grounded metal case
Power:	on/off switch
Power cord:	3-wire UL approved
Weight:	2 lbs.

PLAS ■ LABS, INC. • Lansing, MI U. S. A.

517-372-7177 • fax: 517-372-2857 • www.Plas-Labs.com • Sales@Plas-Labs.com